

# Helium–Neon Laser

A. Simankovich   D. Dedkov

Moscow Institute of Physics and Technology

# Abstract

TODO

# General laser scheme

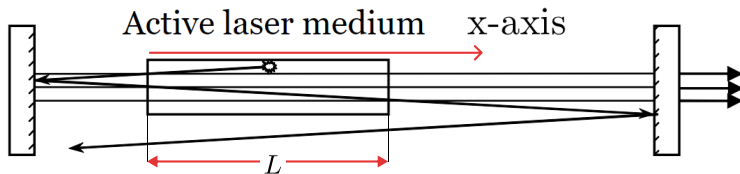


Figure: General laser scheme

# Einstein coefficients


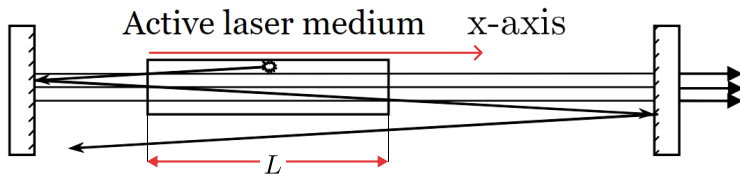


Figure: General laser scheme

# Elementary processes

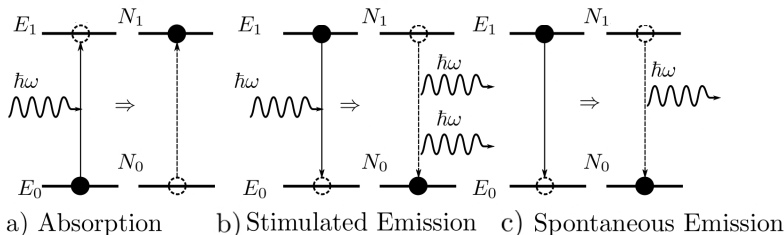

$$\left(\frac{dN_0}{dt}\right)_{\text{abs}} = -B_{01}N_0\rho(\omega) \quad \left(\frac{dN_0}{dt}\right)_{\text{stim}} = B_{10}N_1\rho(\omega) \quad \left(\frac{dN_0}{dt}\right)_{\text{spon}} = -A_{10}N_1$$

Einstein coefficients are the same  $B_{01} = B_{10} = B$

phase, direction and frequency of emitted and external photons are identical.

photons radiate independently in all directions.  $\frac{dN_0}{dt}$  **does not** depend on  $\rho(\omega)$ .

Where  $\rho(\omega)$  – spectral energy density of the isotropic radiation field at the frequency of the transition.

# Gain

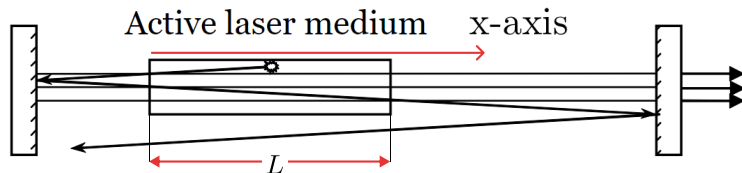


Figure: General laser scheme

Beer–Lambert–Bouguer law states that intensity of light  $I(x)$  changes:

$$I(x) = I_0 \exp(\gamma x),$$

where  $\gamma$  – medium gain coefficient. With length  $L$  gain per period is called **laser gain**:

$$G = \exp(2\gamma L).$$

## Population inversion

The fact that the number of spontaneously emitted photons does not depend on  $\rho(\omega)$  gives us a reason to neglect  $\left(\frac{dN_0}{dt}\right)_{\text{spon}}$  term. Number of photons emitted at a time  $dt$ :

$$\frac{dN}{dt} = \left(\frac{dN_0}{dt}\right)_{\text{abs}} + \left(\frac{dN_0}{dt}\right)_{\text{stim}} = B(N_1 - N_0)\rho(\omega) \quad (1)$$



Therefore:

$$\gamma = \frac{dI}{I} = \frac{dN \cdot \hbar\omega}{\rho(\omega)} = B \frac{\hbar\omega}{v} (N_1 - N_0), \quad (2)$$

where  $v = \frac{c}{n}$  – speed of light inside medium.

$\gamma$  is positive if  $N_1 > N_0$ . This laser principle is called  
**population inversion**

# Polarization of Laser Emission

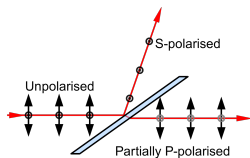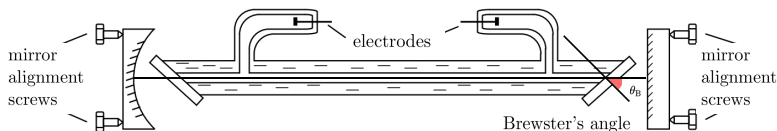


To remove reflections from laser's windows the Brewster's angle properties can be used:

$$r_p = \frac{E_r}{E_i} = \frac{n_2 \cos \theta_i - n_1 \cos \theta_t}{n_2 \cos \theta_i + n_1 \cos \theta_t} \bigg|_{\theta_i = \theta_B} = 0$$



## Malus' law

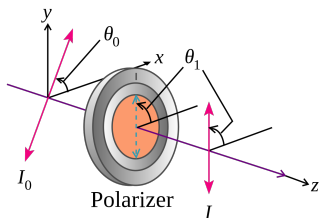

**Figure:** Malus' law (here  $\theta_i = \theta_1 - \theta_0$ )


To study the laser's polarization, polaroid and Malus' law is used:

$$I(\theta_i) = I_0 \cos^2 \theta_i,$$

where  $I_0$  is the initial intensity and  $\theta_i$  is the angle between the light's initial polarization direction and the axis of the polarizer.

# Measurements and Results

# Experimental setup

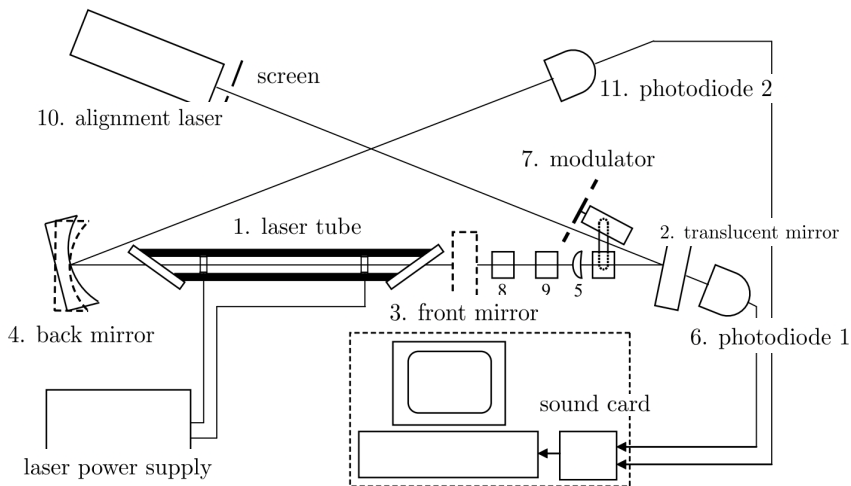


# Experimental Setup



Figure: Photo of laboratory setup

# Laser Gain

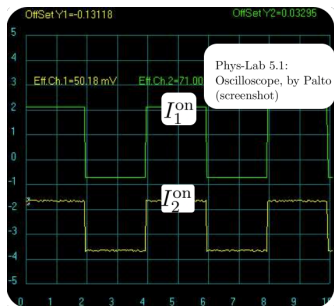

**Figure:** Phys-Lab 5.1:  
Oscilloscope, by Palto  
(screenshot)


Photodiodes are connected to the sound card with the ADC. This gives us a direct way to measure intensity.

$$G = \left( \frac{I_1^{\text{on}}}{I_1^{\text{off}}} \right) / \left( \frac{I_2^{\text{on}}}{I_2^{\text{off}}} \right),$$

where  $I_i^j$  – r.m.s. light intensity.  
Series of measurements gives the following result:

$$G = 1.029 \pm 0.006$$

# Laser Polarization

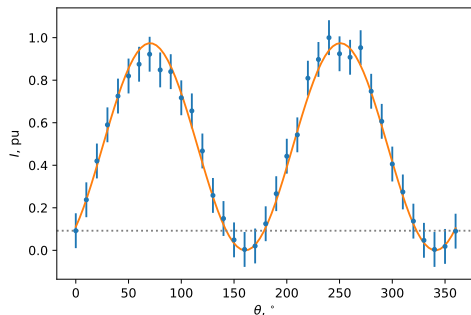

**Figure:** Intensity for different polaroid angles

Interpolating data series in the following form:

$$I(\theta) = A \cos^2 (\Omega\theta + \theta_0),$$

We obtain the following parameters' values:

$$\Omega = 0.998 \pm 0.005$$

$$\theta_0 = (-70 \pm 1)^\circ$$

This demonstrates that Malus's law holds with a great precision.

# Gain

The fact that the number of spontaneously emitted photons does not depend on  $\rho(\omega)$  gives us a reason to neglect  $\left(\frac{dN_0}{dt}\right)_{\text{spon}}$

term. Number of

1. laser tube

laser power supply

4. back mirror

3. front mirror

sound card

2. translucent mirror

6. photodiode 1

11. photodiode 2

7. modulator

screen

10. alignment laser Brewster's angle

mirror  
alignment  
screws

Phys-Lab 5.1:  
Oscilloscope, by Palto  
(screenshot)